

LECKEY



Squiggles+

Perfect blend of style and
function for little ones

Supporting Children

Building on the best of the much-loved Squiggles stander

The incredibly versatile **three-in-one Squiggles stander**, with its minimalist child-focused design, has always been a standout when it comes to happy parents and happy children.

Noted for being **low to the ground**, the Squiggles stander is ideal for encouraging socialization and interaction at eye level, while offering **unsurpassed support, adjustability and easy storage**.

Retaining the best of the Squiggles stander, we have now created the **Squiggles+** with heaps of added benefits so little ones can continue to enjoy standing and be part of the action.



Squiggles+

Perfect blend of style and function for little ones

1

ANATOMICALLY-CORRECT ABDUCTION FOR EARLY HIP JOINT DEVELOPMENT

- Stabilizes the joint and provides optimum loading
- Offering 60° of hip abduction bilaterally in prone, upright or supine
- Promoting proper alignment in the hips, knees and feet as pivot points are in line with the hips
- Hip laterals that shift in sync with the child's leg allowing the hips to abduct and preventing discomfort

2



WRAP-AROUND KNEE SUPPORTS FOR EXTRA COMFORT

- Medial knee supports to distribute pressure, especially suited to abducted standing
- Padded wraps to maintain knee in neutral alignment
- Can be rotated to provide extra lateral support to the knee for low tone children

3

TRANSFER WITH CONFIDENCE

- Larger hip pads to improve carer's and child's confidence and feeling of security when transferring
- Improved pivot base to facilitate full horizontal loading

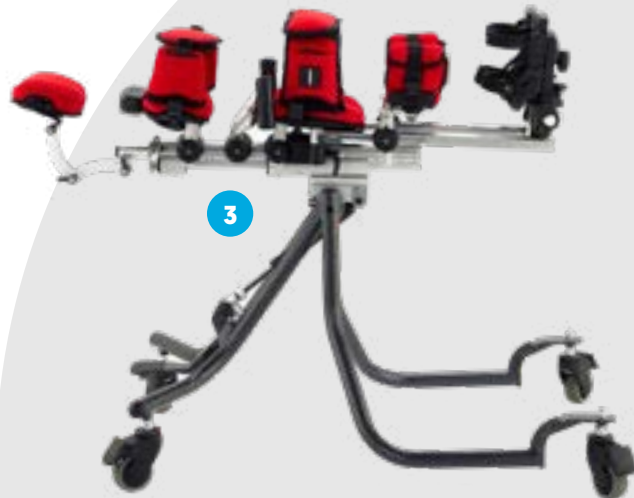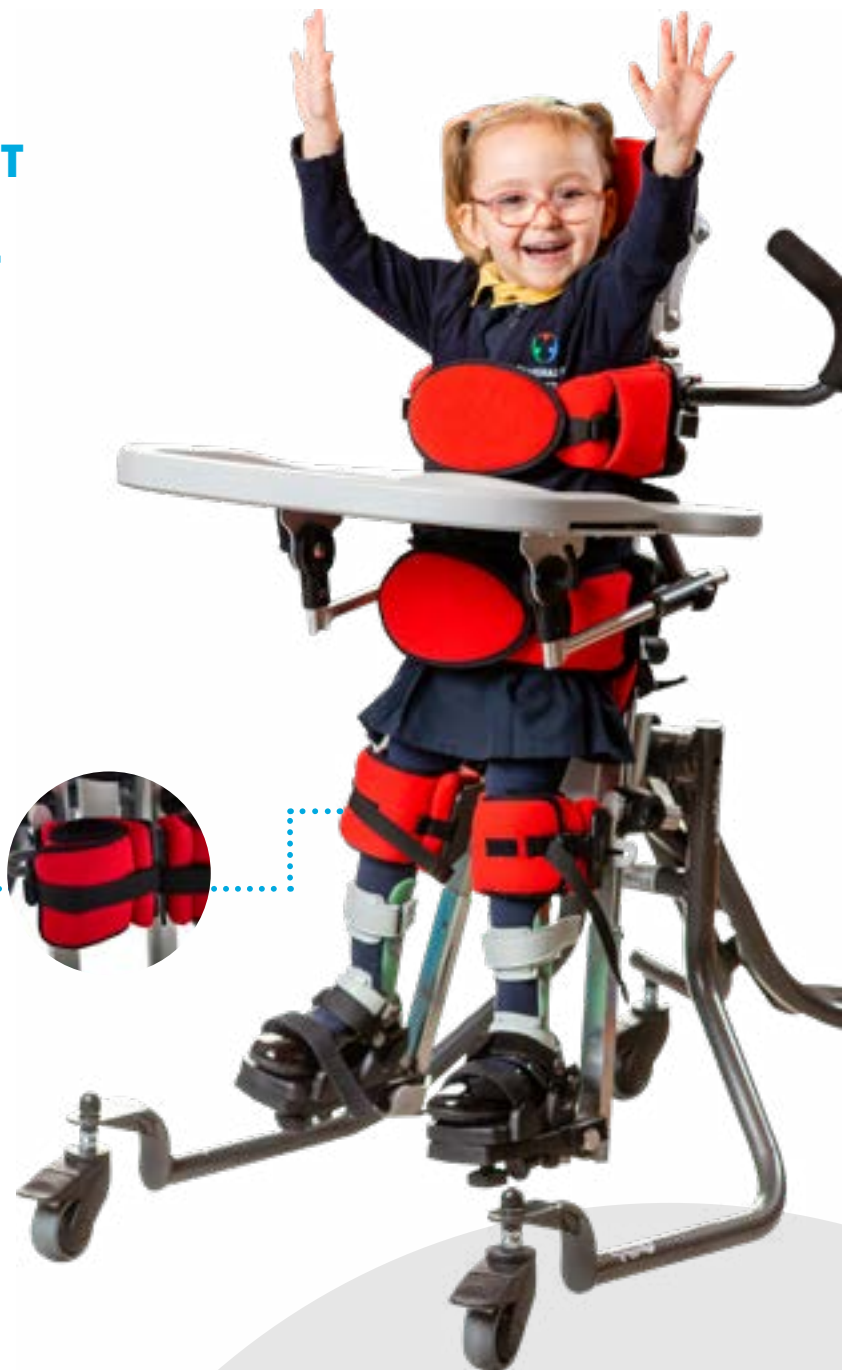

Technical Information

Stander Dimensions

Age (approx.)	1-5
Max User Weight	22kg/48.4lbs
User Height	Min 750mm/28.5 inches Max 1110mm/43.3 inches
Chest Support	Angle adjustment +/- 30° Depth adjustment 50mm/2 inches
Chest Support Height (footplate to armpit)	Min 550mm/21.5 inches Max 820mm/32 inches
Chest Width (distance between laterals)	Min 160mm/6.3 inches Max 230mm/9 inches
Hip Support (height from footplate to mid buttocks)	Min 300mm/11.8 inches Max 530mm/20.8 inches
Distance between midline of Knees	Min 140mm/5.5 inches Max 215mm/8.4 inches
Footplate Angle	Plantarflexion 10° / Dorsiflexion 10°

Tray Size

Tray Height - Prone (footplate to top of tray)	Min 555mm/21.9 inches Max 825mm/32.5 inches
Tray Height - Supine (footplate to top of tray)	Min 480mm/18.5 inches Max 800mm/31.5 inches
Tray Angle Adjustment	Prone 30° / Supine 40°

Support Frame & Pivot Base Dimensions

Stander support frame	Weight 7kg / 16.5 lbs Width 390mm / 15.4 inches Length 680mm / 26.8 inches Height 390mm / 15.4 inches
Pivot Base footprint (unfolded)	Weight 6.5kg / 14.3lbs Width 590mm / 23.2 inches Length 950mm / 37.4 inches Height 515mm / 20.3 inches
Pivot Base footprint (folded)	Weight 6.5kg / 14.3lbs Width 590mm / 23.2 inches Length 660mm / 26.0 inches Height 515mm / 20.3 inches
Pivot Base angle range	Prone 90 - 180° Supine 90-160°